


# Get Started

## DECT indepDant

### Emergency Caller with Assistive Functions

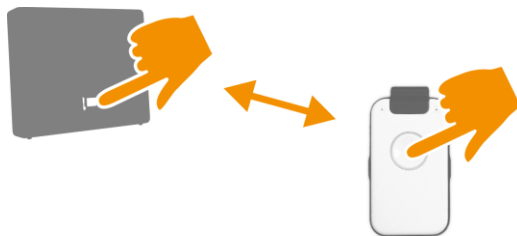

1. Place Device close to DECT base station
2. Activate Registration mode of the base station
3. Press Call Button of Device to start the Registration
4. Pay attention to announcements of Voice Assistant

3


### indepDant overview



**indepDant** is your daily assistant with central Call Button. If you need help, press Call Button and you will be connected with neighbours or friends. Fall Detection or Smoke Alarm automatically trigger Help Call. All activities are conducted by the language assistant. Phone book (10 entries) with list of names is read.

1

### step 3: check telephone connection

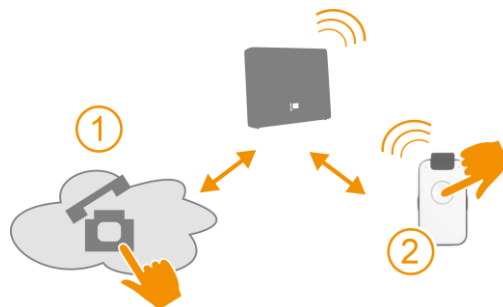

1. Call Device with your mobile phone
2. If Device rings, answer the call


Note: The calling number is stored as the first number to call for help, fall detection and smoke alarm call.

Change or supplement this number with the PC APP.

4

### SET-into-OPERATION

#### step 1: Charge Device

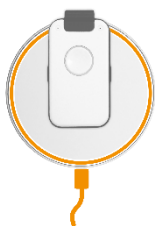

1. Place Charging Disc on a flat, non-slip surface
2. Connect the Charging Disc mains via supplied USB cable
3. Place Device for approx. 30 min on Charging Disc

Note: Device will get a little warm; this is normal

2

### Device is ready to use

Congratulations

now you can

- receive phone calls
- call phone number (see Family & Friends)
- execute a Help Call

following features are enabled:

- Fall Detection with automatic help call
- Smoke alarm call

You can configure via PC APP:

- all contacts for emergency (help, fall, smoke alarm)
- more Family & Friends contacts
- Time and Reminder Messages

5

